


## Rhodamine 6G

## Chemical Properties

CAS No. :	989-38-8
Formula:	C <sub>28</sub> H <sub>31</sub> ClN <sub>2</sub> O <sub>3</sub>
Molecular Weight:	479.01
Storage:	Keep away from direct sunlight Powder: -20°C for 3 years   In solvent: -80°C for 1 year <small>Actual storage temperature shall be subject to the COA.</small>

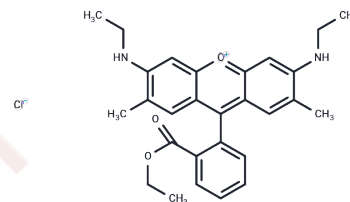

## Biological Description

Description	Rhodamine 6G (Basic Red 1) is a rhodamine analog useful in Pgp efflux assays. It can be widely used in characterizing the kinetics of MRP1- mediated efflux as a laser dye and potential mitochondrial probe.
Targets(IC50)	Others
In vitro	Rhodamine 6G, also known as Rhodamine 590, is widely used as a lasing medium and as a fluorescence tracer. Rhodamine 6G is a fluorescent dye capable of penetrating a living cell. Upon entering the cell, Rhodamine 6G binds to the inner membranes of mitochondria. Based on these observations, it has been proposed that Rhodamine dyes may be used for producing fluorescent images of mitochondria with low background noise and high resolution. Extremely low concentrations of Rhodamine 6G appear to selectively destroy malignant cells in culture, sparing the normal cell populations[2].
In vivo	Compared to their untreated counterparts, Melanoma-transplanted mice receiving Rhodamine 6G demonstrate prolonged survival, inhibited tumor growth, improved clinical parameters and metastases count. It is possible to monitor changes in the interactions between leukocytes and the endothelium by determining the numbers of rolling and adhering leukocytes as well as the total flux of these cells[3]. The Rhodamine-6G enters the circulatory system and labels leukocytes. Twice-a-week 10-6M Rhodamine 6G regimen yield the most prominent results[2].
Cell Research	<ol style="list-style-type: none"> <li>Solution preparation <ol style="list-style-type: none"> <li>Mother solution preparation: Dissolve Rhodamine 6G in an appropriate solvent (such as DMSO or PBS) to obtain the required mother solution, 1-100mM.</li> <li>Working solution preparation: Dilute the mother solution into the required working solution using PBS/DMEM/H<sub>2</sub>O, usually with a concentration range of 1-20µM.</li> </ol> </li> <li>For drug efflux studies <ol style="list-style-type: none"> <li>Cell treatment: Add the dissolved Rhodamine 6G to the cultured cells and incubate with the drug to study the drug efflux process.</li> <li>Efflux kinetics determination: Monitor the fluorescence changes of Rhodamine 6G by fluorescence microscopy or flow cytometry to evaluate the kinetics of drug efflux.</li> </ol> </li> <li>As a mitochondrial probe: <ol style="list-style-type: none"> <li>Dye working solution concentration: Usually the working solution concentration ranges from 1-5 µM.</li> <li>Cell incubation: Add the dissolved dye to the cell culture medium and incubate for a</li> </ol> </li> </ol>

Cell Research	<p>certain period of time (usually 30 minutes to 1 hour) to allow it to enter the cells and accumulate in the mitochondria.</p> <p>3. Fluorescence microscopy observation: Use a fluorescence microscope to detect the fluorescence signal of mitochondria and observe the distribution and function of mitochondria.</p> <p>3. As laser dye:</p> <p>1. Dye working solution concentration: Usually the working solution concentration is 10-100 <math>\mu\text{M}</math>.</p> <p>2. Laser dye application: Apply the dye to the laser system, or use it for fluorescence microscopy or flow cytometric analysis.</p> <p>3. Fluorescence signal detection: Detection is performed by fluorescence microscope, flow cytometer or other laser analysis equipment.</p> <p>The above information is based on published literature. Experimental procedures should be appropriately modified to meet specific research demands.</p>
---------------	--

### Solubility Information

Solubility	<p>H<sub>2</sub>O: 83.3 mg/mL (173.9 mM),Sonication is recommended.</p> <p>DMSO: 15.625 mg/mL (32.62 mM),Sonication is recommended. (<math>&lt; 1</math> mg/ml refers to the product slightly soluble or insoluble)</p>
In vivo Formulation	<p>10% DMSO+90% Saline: 1.56 mg/mL (3.26 mM),Suspension.</p> <p><i>Please add the solvents sequentially, clarifying the solution as much as possible before adding the next one. Dissolve by heating and/or sonication if necessary. Working solution is recommended to be prepared and used immediately. The formulation provided above is for reference purposes only. In vivo formulations may vary and should be modified based on specific experimental conditions.</i></p>

### Preparing Stock Solutions

	1mg	5mg	10mg
1 mM	2.0876 mL	10.4382 mL	20.8764 mL
5 mM	0.4175 mL	2.0876 mL	4.1753 mL
10 mM	0.2088 mL	1.0438 mL	2.0876 mL
50 mM	0.0418 mL	0.2088 mL	0.4175 mL

Please select the appropriate solvent to prepare the stock solution, according to the solubility of the product in different solvents. Please use it as soon as possible.

Note: The dilution table applies only to solid products. For liquid products, please calculate the stock solution based on the stated concentration and/or density.

### Reference

- Wang J,et al. Uniformly aligned Ag NPs/graphene paper for enhanced SERS detection of pesticide residue. Spectrochim Acta A Mol Biomol Spectrosc. 2025 Mar 15;329:125535.
- Liu Y, Zhou M, Sun J, et al. Programmed BRD9 Degradation and Hedgehog Signaling Activation via Silk-Based Core-Shell Microneedles Promote Diabetic Wound Healing. Advanced Science.2024: 2404130.
- Yang J,et bal. Monitoring molecule translocation through plasmonic nanopores based on surface enhanced Raman scattering. Spectrochim Acta A Mol Biomol Spectrosc. 2025 Mar 15;329:125511.
- Zhou R,et al. Magnetic metal-organic frameworks-based ratiometric SERS aptasensor for sensitive detection of patulin in apples. Food Chem. 2025 Feb 28;466:142200.

**Inhibitor · Natural Compounds · Compound Libraries · Recombinant Proteins**

This product is for Research Use Only · Not for Human or Veterinary or Therapeutic Use

Tel:781-999-4286 E\_mail:info@targetmol.com Address:34 Washington Street,Wellesley Hills,MA 02481